

# NEPLUS (NE-SM1) SUPER MESH NEBULIZER

- Vibrating mesh technology (Nickel Free)
- Rechargeable battery with smart phone adapter  
(No battery change for more than 5 years)
- Noiseless operation
- Preventing infection by boiling

**NEW**



model **NEPLUS (NE-SM1)**

## Features

- Ultrasonic type nebulizer that sprays liquid medication into lung as aerosols.
- Ultrasonic vibrator mesh releases medication via 5 $\mu$ m of fine holes.
- Very small size and noiseless Mesh technology.
- Very convenient to use and to carry.
- Inhalation is possible at over 75° of tilt angle.
- Long operating time under rechargeable battery.


CE  
0120



# NEPLUS (NE-SM1) SUPER MESH NEBULIZER



## Technical Data

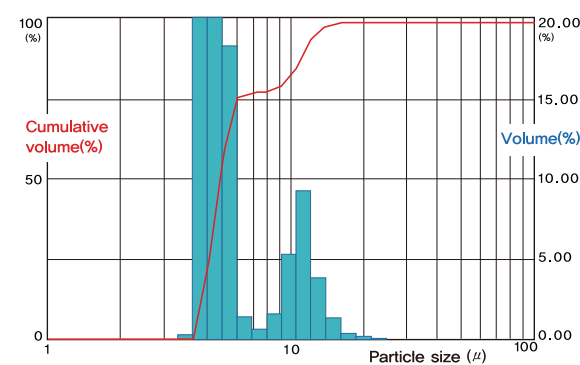

Type	Pulmicort	Ventolin
Nebulizer Type		
Ultrasonic Nebulizer	Not Good	Good
Mesh Nebulizer	Good	Good

## Specifications


Model	NE-SM1
Power & Consumption	Lithium ion battery 3.6~4.2V DC 1.3 W
Dimensions	(W)×(D)×(H) = 52mm×53mm×109mm
Weight	Main unit : 88g, Mask : 50g
Nebulization rate	Approximately 0.25ml/min
Particle size	Average 5 $\mu$ m
Volume of medication bottle	6ml
Operating time under full charge	5 hours 30 min
Charging life cycle of the battery	over 500 cycles
Ultrasound frequency	Approximately 108 KHz
Operating environment	10~40°C , related humidity 30~85%
Battery charging environment	0~45°C
Storage environment	-20~60°C
Contents	Nebulizer unit, Adult mask(optional), Child mask, Mouthpiece, Instruction manual, Carrying bag

\*All specifications are subject to change prior notice.



KTMED INC. #610, Hanshin IT Tower, 272 Digital-Ro Guro-Gu, Seoul, Korea.  
Tel : +82-2-2108-7890 Fax : +82-2-2108-7894  
E-mail : master@ktmed.com

## Composition

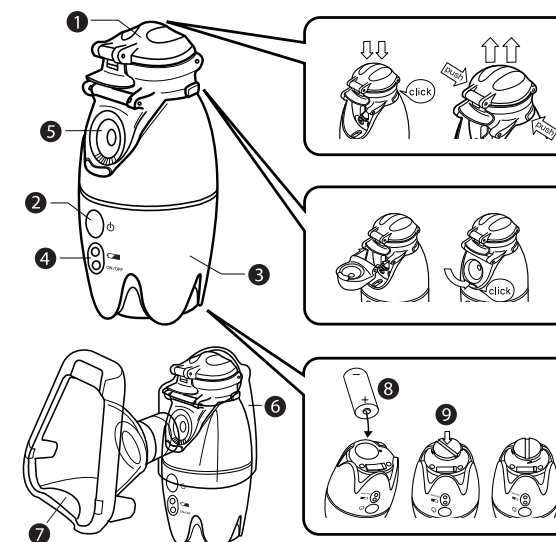

## Components

Main unit	Mesh cartridge	Medication bottle
Mask adapter	Small inhalation mask	Big inhalation mask(Optional)
Rechargeable lithium battery	Mouthpiece	AC power cable (Individual purchase)
		Smart phone adapter

\*For repairing or repurchasing of parts, please call our company or local distributor.

